

AmaMedic

AM-4610Eye Massager

User Manual



Relax Music



Graphene Heating



Air Compression



Vibrating




180° Foldable

This user manual provides both operation information of this product and its detailed specifications. Read this manual and understand it thoroughly before operation.

Operation Guide


Bluetooth Version

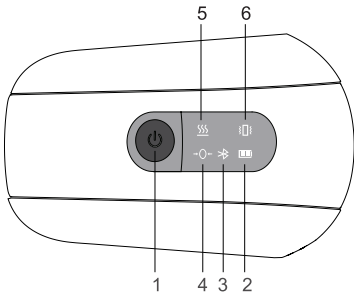

- 1, Power : Long Press ON/OFF
Short Press Switches Modes
- 2, Battery: LED Shining with no power
LED Breathing when Charging
LED Lighting when full charging
- 3, Bluetooth
- 4, Air Compression
- 5, Heating
- 6, Vibrating

Mode	Air compression	Heating	Vibrating	Bluetooth	Timer
Energy Mode	√	√	√	√	10
Relax Mode	√	√		√	15
Relieve Mode	√	√		√	15
Sleep Mode	√			√	20

Music Version

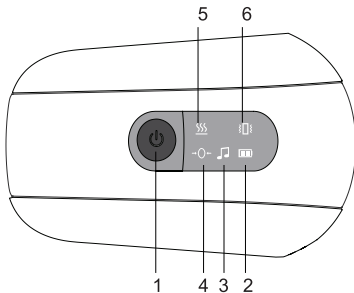

- | | |
|---|---|
| <p>1, Power : Long Press ON/OFF
Short Press Switches Modes</p> <p>2, Battery: Shining with no power
Breathing when Charging
Lighting when full charging</p> | <p>3, Music</p> <p>4, Air Compression</p> <p>5, Heating</p> <p>6, Vibrating</p> |
|---|---|

Mode	Air compression	Heating	Vibrating	Music	Timer
Energy Mode	√	√	√	√	10
Relax Mode	√	√		√	15
Relieve Mode	√	√			15
Sleep Mode	√			√	20

Product Function

- 1, Buzzing Prompt
- 2, Multiple Vibrating Function
- 3, 104°F Heating
- 4, Single button operation
- 5, USB Type-C jack
- 6, Power shortage light flash, Buzzer sound "DI" per 2seconds, When no power it shut down, when power on it "DI" 3times alarm.
- 7, Four Modes: Energy Mode, Relax Mode, Relieve Mode, Sleep Mode

If you

- 1, Are pregnant
- 2, Are shortsighted over than 800 degrees
- 3, Are seriously glaucoma, cataract
- 4, Are under 7 year-old
- 5, Are sequela of apoplexy
- 6, Are eye injury, in the period of rehabilitation of finished surgery (double eyelid)
- 7, Were Glaucoma, cataract, retinal detachment and other serious ocular region surgery
- 8, Are suffering from serious eye inflammation

We advise you to consult with your doctor before using this device.

Safety Note

Attention: This device is designed for private use, not intended for medical or commercial purpose. Children should be supervised to ensure that they do not use this device as a toy.

- 1, Suggest to Adjust elastic suitable and Sit leaning against the sofa or chair when using. That will be more comfortable to use.
- 2, Suggest to use this device one or two times per day, Each time should be 10-20 minutes. And do not open your eyes when using this device.
- 3, Please using this device after charging.
- 4, Please cleaning surface of this device with wet towel or Cotton ball with 75% medical alcohol.
- 5, Do not put this device close fire,corrosive liquid and other adverse environments.

Product Specifications

- 1, Name: AM-4610 Eye Massager
 - 2, Input: DC5V/1A
 - 3, Rated Power: 5W
 - 4, Battery : Built-in Li-ion battery
3.7V 1200MAH
-

AmaMedic