# PhysioPedal<sup>®</sup>

User Manual EN ES FR 日本 中文

by Nobol<sup>®</sup>

*Plug-In / Resistance*

**Hello**

**Hola**

**Bonjour**

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# Contents

Safety instructions .....	5
What's in the box .....	6
Specifications .....	7
Features .....	8
Remote Control / Display .....	9
Remote Control Buttons .....	10
Getting Started .....	11
Resistance Mode .....	12
Assisted Mode .....	13
Troubleshooting .....	14
Warranty .....	15

# Welcome

# Safety Instructions

When using The PhysioPedal, basic precautions should always be followed, including the following:

**Read all instructions before using The PhysioPedal.**

## Warning

To Reduce the risk of burns, fire, electric shock, or injury to persons:

Disconnect the charger from the PhysioPedal when not in use, and before putting on or taking off parts.

Use The PhysioPedal only for its intended use as described in this manual.

Do not use attachments not recommended by the manufacturer.

- Never operate The PhysioPedal if it's damaged, or not working properly, if it has been dropped or damaged, or dropped into water.
- Return The PhysioPedal to a service center for examination and repair.
- Never drop or insert any object into any opening.
- Do not use it in the rain.
- Do not operate where aerosol (spray) products are used or where oxygen is administered.
- For household use only.

This appliance can be used by children aged 8 years and above and by individuals with reduced physical, sensory, or mental capabilities, or a lack of experience and knowledge, if they have been given supervision or instruction regarding the safe use of the appliance and understand the hazards involved. Children should not play with the appliance. Cleaning and user maintenance should not be performed by children without supervision. The appliance is intended to be used only with the power supply unit provided with it.

# What's in the box

**A** PhysioPedal Body

**B** Base

**C** Wireless Remote

**D** Plug-In Cable & Adapter

**E** Oversized Pedals x 2

**F** Handles x 2

**G** Base Screws x 4

**H** Wrench


# Specifications

<b>Product Model:</b> PP0002	<b>Net Weight:</b> 8.61 Kg.   19 Lbs.
<b>Rated Power:</b> 72W	<b>Gross Weight:</b> 11.25 Kg.   24.8 Lbs.
<b>Product Size:</b> 20.2 x 10.6 x 14.6 in.	<b>Body:</b> 5.4 Kg.   12 Lbs.
<b>Output:</b> 24V=3A	<b>Handles:</b> 0.3 Kg.   0.7 Lbs.
<b>Input:</b> 110~220V AC	<b>Comfort Pedals:</b> 1.6 Kg.   3.6 Lbs.

## Cleaning and Storing


Make sure the power is **OFF**. Use a damp, clean cloth and gently wipe your PhysioPedal. Store in a clean, cool, dry place, not in direct sunlight, when not in use. Any other servicing should be performed by an authorized service representative.

# Features

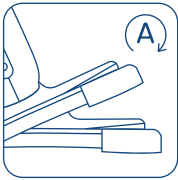
## Tilt-sensitive safety feature

10° angle triggers alarm and auto-stops motor
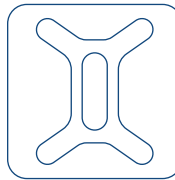
## Ergonomic Base

Enhanced stability, lightweight design

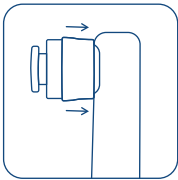
## Quick Release Pedal

Pull the shaft and easily insert or release the pedal

## Remote Control

Senior friendly Intuitive Interface
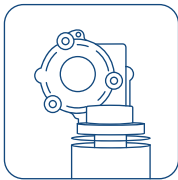
## Motor

Medical grade DC electric motor


## Interactive Screen

Tracking and reporting of all usage, calories, distance, time used

# Remote Control / Display

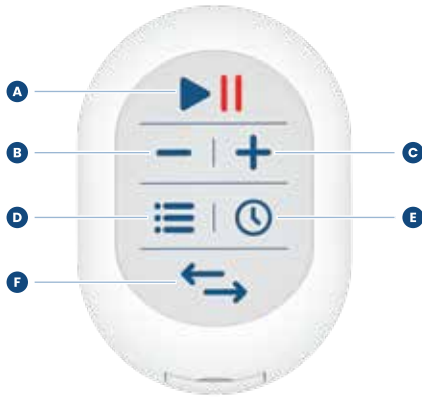- A** Start/Stop
- B** Decrease Speed
- C** Increase Speed
- D** Menu
- E** Timer
- F** Switch Direction

- A** Rotations per minute
- B** Timer
- C** Timer Countdown
- D** Calories
- E** Steps

# Remote Control Buttons

Press the **Start/Stop** button once to start the Pedals or Handles to move.

Press **Start/Stop** button again to stop.

The unit will start a slow speed of 30 RPM, to increase the speed press the **Plus** button.

To decrease the speed press the **Minus** button.

To set the timer, press the **Clock** button.

The PhysioPedal will start at 10 minutes.

Each press will increase the time by 10 minutes up to 90 minutes total.

To change between Time, Calories, and Steps, press the **Menu** button.

To switch the rotational direction of pedals, push the **Two Arrows** button.

---

If the remote control stops working, you may need to reconnect the remote control to your PhysioPedal. Turn the PhysioPedal **ON**  then, on the remote control, press and hold the **Switch Direction** button  for 3 seconds until you hear a beep.

# Getting Started


## Assembly

To assemble the PhysioPedal, use the **Tool** (pic1) to screw the **4 Screws** to attach **the base to the body** of the PhysioPedal.

Attach **the Pedals or Handles** into the Cranks. Switch from Pedals to Handles with a simple action, simply pull the knob back and slide the Pedals or Handles in or out (pic2).

Plug the **Power Cord** (pic3) into the Front of the Unit, screw until hand tight.


Put the **Remote Control** (pic4) over the magnetic screen.

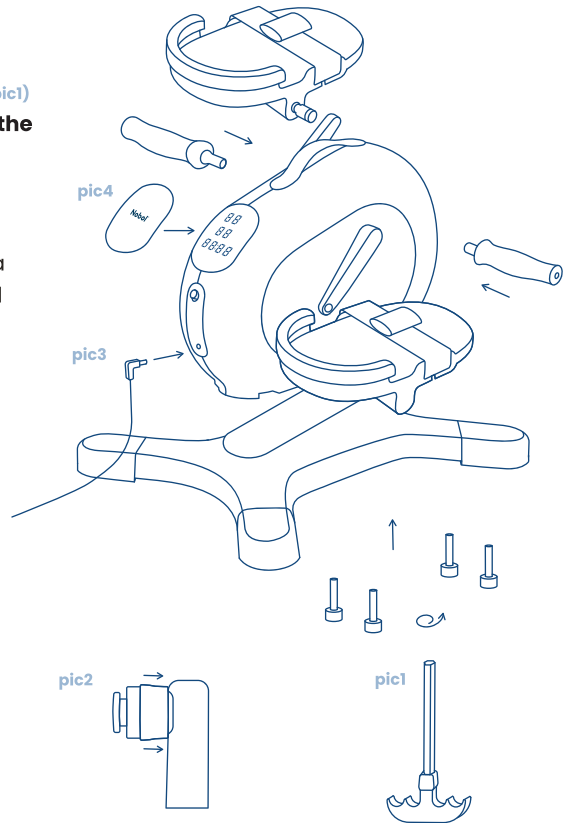
# Resistance Mode

(Default Mode)


To **activate** the **PhysioPedal**, ensure it is connected to an **electrical outlet**. Press the **Power Button** located on the front of the unit to initiate **power**.

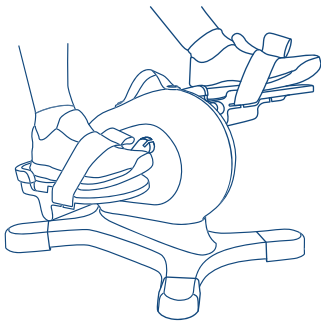To **increase resistance**, simply use the **Remote Control's + button** and select your preferred **resistance level**.


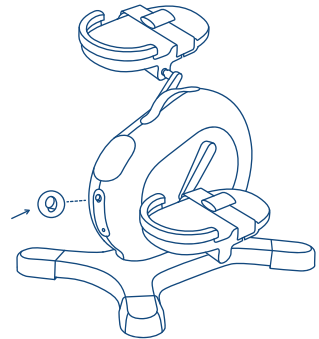
When activating the **PhysioPedal**, start at **Level I**. **Sit comfortably** and position both feet on the **pedals**. **Begin pedaling**.
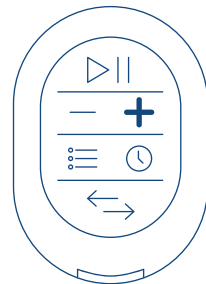
# Assisted Mode

To activate **Assisted Mode**, simply press the **play button** on the **Remote Control**.

To increase speed in the assistance level, press the **+** button on the **Remote Control** to your preferred level. **Sit comfortably** and position both feet on the **pedals**. **Begin pedaling**.

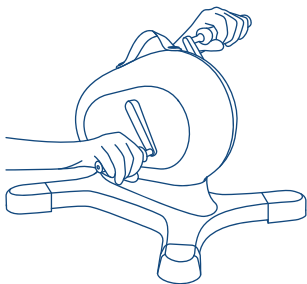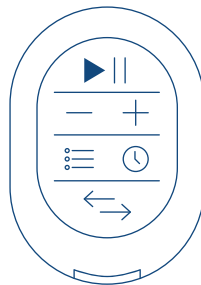For upper body exercise, place the **PhysioPedal** on an elevated platform or table, and position your hands on the grips.

Upon completion, simply turn off the PhysioPedal by pressing the **Switch OFF** button.

# Troubleshooting

If you encounter any issues while using this product, please consult the provided tips and troubleshooting methods. If the problem persists, power off the device

<b>Use Failure</b>	<b>The reason</b>	<b>Solution</b>
Product Not Functioning (Cannot be Turned On)	The power switch is in the off state	Put the power switch on
	The adapter is not powered on	Power on the adapter and check if the DC interface is in good contact
Remote control not working	The battery insulation sheet has not been removed	Remove the insulation sheet
	The battery of the remote control is depleted	Replace with a new battery
After turning on, the buzzer will sound an alarm and display	The product is tilted or has not been positioned on a flat surface	Place the product on the flat ground
Product Stops during operation	Incorrect setting for the timer	Press the power button to restart
Stop working during use, alarm display	When the pedaling speed exceeds the set RPM, an overspeed alarm is activated, which is part of the normal overspeed alarm function	Reducing pedal speed will restore normal operation

# Warranty

**contact our customer service by emailing us at [cs@Nobol.com](mailto:cs@Nobol.com) to start your warranty period.**

If within one year from the date of purchase, this device fails due to a defect in material or workmanship, Nobol Inc. will repair or replace the device or necessary components, free of charge.

## **This warranty excludes:**

- A. Damage caused by accident, abuse, mishandling, or transport
- B. Units subjected to unauthorized repair
- C. Damage exceeding the cost of the product
- D. Deterioration of the delivered product resulting from abnormal storage and/or safeguarding conditions on the client's premises
- E. Products taken apart or repaired by anyone other than an authorized service representative
- F. Failure to provide the dated proof of purchase and serial number on your device
- G. Products purchased from an unauthorized retailer
- H. Please do not try and repair the unit yourself, as this will void your warranty

\*Some states and countries do not allow a limitation of damages, so the foregoing limitation may not apply to you. This warranty guarantees specific legal rights, but other rights may vary between countries and from state to state. This warranty service is available through the operating policies and procedures of Nobol Inc., and all Nobol international distributors. If you are a client, please contact the appropriate Nobol business entity dealer, or retailer, from whom you directly purchased the device for warranty and return authorization procedures. This warranty is initiated and executed by Nobol Inc.

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Designed with ♥ in California by **Nobol, Inc.** Assembled in China