

WARRANTY CARD

Model	Relax FIR Sauna
Model No.	FIR9008E-2sw
Serial No.	G&MSA09311
Name	
Address	
Telephone	
E-mail	
Expiration Date	<u>One year from date of purchase</u>
Purchase Date	
Distributor Stamp	
<p>1. We sincerely support the best service to you. Please complete the warranty card and retain with operating instruction.</p> <p>2. Under our terms and conditions we cannot provide free repair or service under the following circumstances:</p> <ul style="list-style-type: none">❶ Damaged caused by improper use (i.e. not following instruction manual.)❷ Damaged caused by the user changing or repairing the product without prior consent from the manufacturer.❸ Damaged caused by accident, fire, flood or any other natural disasters.❹ <u>Warranty is void if proof of purchase is missing.</u> <p>3. Please keep the warranty card safe, as it cannot be re-issued.</p>	

Relax FIR Sauna



Operating Instruction

- ❖ **Elderly People, Children and those who require medical care must use Relax Sauna under supervision of another adult!**
- ❖ **Elderly People, Children and those who require medical care must consult with physician before using the Relax Sauna!**

- Please read the instruction book carefully before using the appliance.
- The instruction book should be retained with the appliance for future reference.




CAUTION:


- Beware of the following safety precautions on the electric cord to avoid accident:
 - ① Use AC 120V 15A outlet and do not share the outlet with other electrical equipment.
 - ② When the appliance is in operation, do not twist, pull, or bend the electric cord. Do not put any heavy object on the cord.
 - ③ Once the cord is torn or faulty, do not use the appliance and send to the company for repair.
- For cleaning and maintenance:
 - ① Use soft and damp cloth or mild detergent to wipe. Do not use volatile cleaning liquid.
 - ② After a period of usage, use a powerful hair dryer (use only cool wind) to spout out the dirt in the FIR radiators.
 - ③ The FIR radiators' floor mat can be removed for washing.
- It is dangerous to make any alterations or modifications to the appliance's specifications.
- If a fault occurs please turn off the switch on the controller, and then remove the plug from your power supply. Contact your supplier immediately for advice.

STORAGE


		<p>1. Take out the chair, then shrink the frame and take out the heater. (Before you take out heater shrink the frame first, then take it out and take off the tubes – see diagram on left)</p>
		<p>Note: When shrinking frame use right hand to press the “lock button” and use your left hand to slide the heater inward (see diagrams on left).</p>
<p>2. Fold cover as show (see diagrams above): Pull in the right-side → Fold the bottom in → Pull left-side in → OK.</p>		
		<p>3. Complete folding cover as shown: with hoop & loop tape face up, turn point A up and twist cover into “8” shape. (Clean cover internally with damp cloth, dry before folding cover)</p>
<p>4. Bend from point B to point C. → Fasten velcro strap → Completed.</p>		

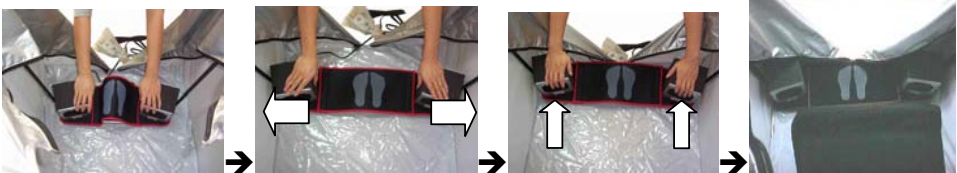


6. Pull the two radiators apart until you hear a “click” sound. Push the radiators near to the zipper. Place chair inside.

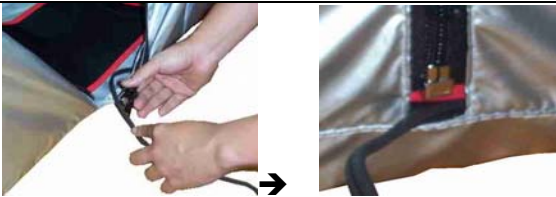

7. Pull the cable out under the zip → Attach the zip → Zip up.



8. Adjust the right tube-tip right below seam and fix the base of controller on the seam and tip of tube together.



9. You are now ready to use your Relax Sauna. Please read the instructions below.

【HOW TO USE】

1. Read warnings on page 3.
2. Install the sauna as shown above then switch on the controller and set the timer.
3. The user goes into the sauna. Zip up the cover
4. Rest feet on the foot mat. 【Fig. 1】
5. Drink plenty of water before and after using Sauna.
6. Suggested usage duration: About 15 min
7. For Personal hygiene, place a towel around the neck. 【Fig. 2】



【Fig. 1】




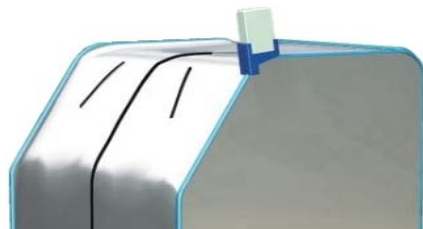



【Fig. 2】

SPECIFICATIONS

Model:	CH-9008	Size:	80x80x115 cm
Power source Voltage :	<input type="checkbox"/> AC 100V, 60 Hz <input type="checkbox"/> AC 110V, 60 Hz <input type="checkbox"/> AC 120V, 60 Hz <input type="checkbox"/> AC 220V, 50 Hz <input type="checkbox"/> AC 230V, 50 Hz <input type="checkbox"/> AC 240V, 50 Hz	Power:	1500 W +5%-10%
Cord	Plug to controller: 232 cm Controller to heater: 245 cm	Weight:	Cover: 2 kg Heater: 6 kg Chair: 3.2kg

PARTS IDENTIFICATION

 <p>Reflect Cover</p>	 <p>FIR heater and tubes</p>	 <p>Chair Max. Load: 350 lb</p>
 <p>The base of controller is set up as shown</p>		 <p>Controller Consists of a Timer and 2 Switches. Each radiator has it's own switch.</p>

IMPORTANT SAFETY INSTRUCTIONS

- Use the appliance on an even floor and beware of sharp-pointed objects that can damage the FIR cover.
- Remove all metal objects (Necklace, watch, ring...etc) before use.
- As Relax Sauna is an electric appliance, **Do not** use it in places with high humidity or direct water contact (such as a bathroom) or touch controller with wet hands, so as to avoid electric shock.
- For personal hygiene, spread a towel on the chair and tie a small towel round the neck position when use.
- Use the Main-Switch to turn on and off the power. The timer is to turn clockwise to the time desire. **Do not** turn counterclockwise. Let the timer stop automatically.
- **Do not** touch the face of the radiator and the metal part of the chair, as they are hot.
- Avoid Tripping, please zip down the reflect cover completely after use.
- If feel the heel is too hot, please move the Sauna away from the wall.
- Usage time is based on individual comfort and to be adjusted accordingly. The recommended time is 10~15minutes per session.
- **Do not** use the Relax FIR Sauna under the following conditions:
 *Pregnant women *Infant *Under the influence of alcohol
 *Under a doctor's instructions not to use sauna
 *Suffering from deep cut, or Fresh wound.
- For the elderly and chronic disease patients, please use the sauna under supervision of an adult and keep the session short (4~6 minutes) until the user is accustomed to the Relax Sauna.
- People with special medical problems should consult doctor before use.
- At any point, if you feel very uncomfortable, stop using immediately, rest, and consult doctor if necessary.

INSTALLATION

1. Place cover on floor, release Velcro to release the twisted cover.

2. Open the cover fully, then pull the zip to the bottom and unfasten it, as shown in the diagrams on the left.

3. Take out the FIR heater and the tubing. Set up the tubes and fix the tubes onto the heater as shown in the *diagrams above*. Please note that the left tube and right tube are different – *see diagram above*.

4. Put the Radiators inside the cover as shown in the *diagrams on the left*.

5. Put the tubes inside the cover as shown in the *diagram on the left*.