

# JOBST® Bella™ Strong

Ready-to-Wear Armsleeve and Gauntlet



**New dimensions  
in effectiveness  
and comfort**

True confidence  
for patients with  
lymphedema and  
edema



#1 physician recommended



Comfort, Health and Style!

# Improved quality of life

Lymphedema can be one of the most troubling complications after breast cancer treatment. Optimal fitting compression garments effectively control lymphedema and edema in the upper limb, giving true confidence and enhancing quality of life.

## Improved **JOBST® Bella™ Strong** – Uncompromised medical efficacy for every patient

For over 50 years BSN has developed technologies addressing patients needs. New, improved **JOBST® Bella™ Strong** is the result of our continual patient-focused research. It offers an unrivaled range of patients with highly effective compression therapy with all-day wearing comfort.

**JOBST® Bella™ Strong** Ready-to-Wear Armsleeve and Gauntlet – worn together or separately – effectively control moderate to severe lymphedema and edema.

## Improved comfort Better compliance Enhanced quality of life

Patients are more likely to be compliant when they feel comfortable wearing a compression garment. The advanced **JOBST® Bella™ Strong** two-way stretch armsleeve has been developed to make compression therapy more comfortable without compromising medical efficacy.

## Comfort and quality

### **JOBST® Bella™ Strong – Armsleeve**

- **Comfort Flat wristband** reduces rolling and improves durability
- **Improved elastic silicone band\*** with a true 2" footprint
  - Provides secure hold up
  - Helps prevent curling edges
  - Helps prevent constrictions
  - Color matched to gauntlet

\* Available with knitted top band.

### **JOBST® Bella™ Strong – Gauntlet**

- **Significantly flattened seams** avoid thumb web
- **Extended wrist** avoids pinching at the wrist
- **Comfort Flat wrist seam** reduces rolling




# JOBST® Bella™ Strong

## Naturally effective

The new **JOBST® Bella™ Strong** provides patients with the latest development in fiber science. The skin-friendly multifiber nylon combines the clinically effective management of lymphedema with a natural feeling. All day long, day after day.

**Bella™ Strong** features excellent wicking properties with a cotton-like feel. **JOBST® Advanced comfort™** helps to regulate relative humidity, feeling cooler in warm conditions and warmer in cool conditions.

### JOBST® Advanced comfort

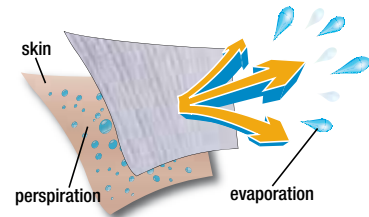

- Quickly wicks moisture
- Keeps you cool and comfortable
- Feels soft to the skin

### Superior Comfort Superior Performance

Due to its outstanding moisture management properties it takes Bella™ Strong only seconds to wick away moisture.\*

**In addition, the open knit structure allows for an improved air permeability that keeps patients' skin drier than the competition.**



Bella™ Strong: 20-30 mmHg



Bella™ Strong: 30-40 mmHg



Competitor I: 30-40 mmHg



Competitor II: 30-40 mmHg

All pictures were taken at equal (50x) magnification

\* Laboratory testing results are intended to illustrate a product's performance under controlled conditions. Actual use results may vary.

“

The sleeve is so much more comfortable than my heavy one. It is much easier to put on and stays in place.

”





# Engineered for living

## Better Then Ever

Helping you manage your lymphedema is a top priority and the major focus here at JOBST®. The fabric in JOBST® Bella™ Strong Armsleeves has been newly engineered for the daily challenges and demands of an active lifestyle.

## Competitive Testing: Puncture Resistance

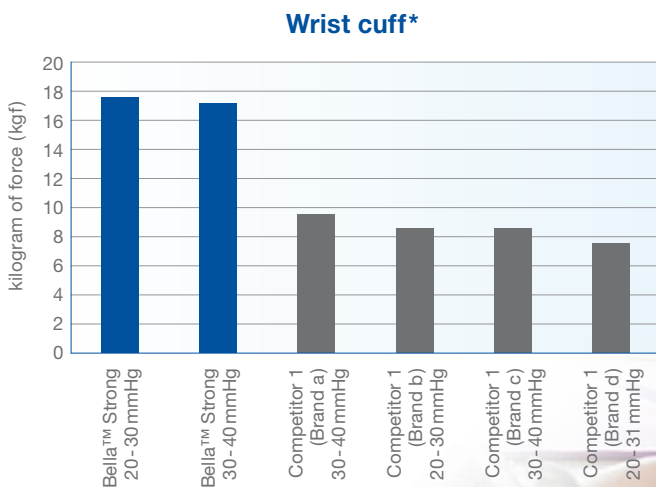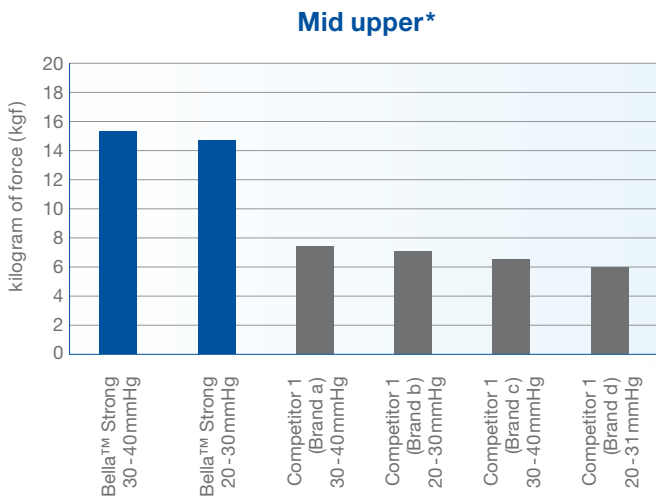

\* Laboratory testing results are intended to illustrate a product's performance under controlled conditions. Actual use results may vary.

“ The armsleeve doesn't stretch out after wearing for a few days. The compression on arm is more even, so no area becomes more swollen than another. ”



## Innovative sizing system – Covers an unrivalled range of patient dimensions

Considering over 50 years of scientific research and practical experience as well as the data from over 36,000 lymphedema patients for the armsleeve and over 33,000 for the gauntlet, JOBST® has developed an innovative measuring-size-system.




Catering to all relevant armsleeve measuring points – wrist, mid lower and axilla – this thorough approach, accepts the largest range of patient dimensions. **JOBST® Bella™ Strong** armsleeve sizing begins at the forearm, as data shows that we can accommodate more patients; this allows for multiple size options in the wrist and axilla for greater coverage, recognizing that every patient is unique.


### JOBST® Bella™ Strong Armsleeve: Ten-size system provides over 80 % fit coverage

Competitors' wrist centric size determination limits fit coverage.

Fit of 36,994 patient measurements

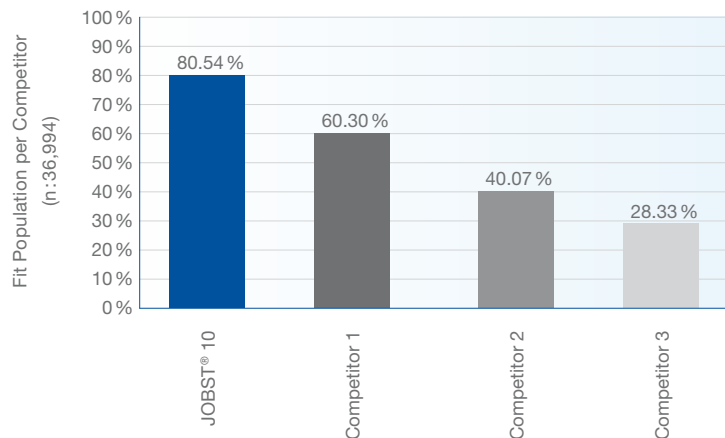

Fit of 33,748 patient measurements

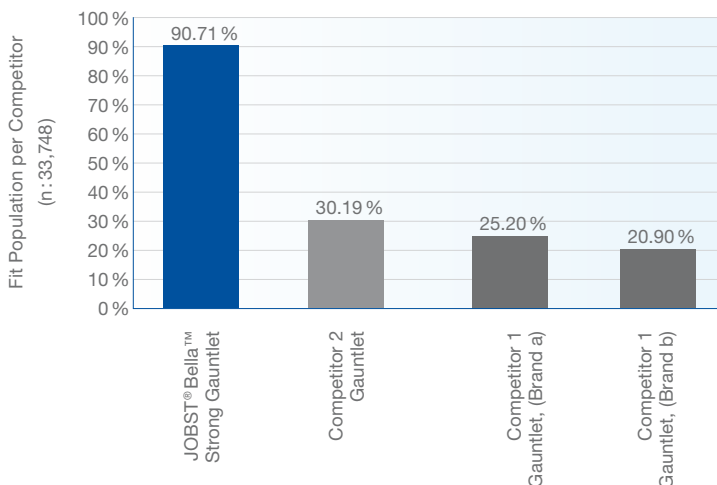

### JOBST® Bella™ Strong Gauntlet: Over 90 % fit coverage

With the new eight-size system the **JOBST® Bella™ Strong** gauntlet reaches an outstanding range of fit coverage.

# Improved for your needs

## Sizing: Armsleeve

**JOBST® Bella™ Strong** Armsleeves are designed with overlapping sizes that include measurements at the wrist, mid forearm and axilla so patients get an optimal fit without compromising medical efficacy.

Size	Wrist (C)	Mid forearm (D)	Axilla (G)
1	14 - 18 cm	17 - 21 cm	22 - 31 cm
2	14 - 18 cm	17 - 21 cm	31 - 41 cm
3	14 - 18 cm	21 - 25 cm	25 - 34 cm
4	16 - 19 cm	21 - 25 cm	28 - 38 cm
5	17 - 20 cm	21 - 25 cm	34 - 43 cm
6	14 - 18 cm	25 - 29 cm	29 - 38 cm
7	17 - 20 cm	25 - 29 cm	33 - 43 cm
8	18 - 21 cm	25 - 29 cm	38 - 48 cm
9	16 - 19 cm	29 - 33 cm	32 - 42 cm
10	19 - 22 cm	29 - 33 cm	36 - 46 cm

### Even more improvement:

All sizes are available in 2 lengths:

**Regular** (up to 43 cm)

**Long** (greater than 43 cm)



## Sizing instructions

Finding the right size for each patient is easy since stepwise exclusion leads to the perfect solution:

1. Measure mid lower circumference. This leaves 2 or 3 choices.
2. Measure axilla circumference. This leaves 1 or 2 choices.
3. Measure wrist circumference. This leaves 1 or 2 choices.

If there are still 2 options left, we recommend choosing the larger size for ease of donning, overall comfort and ultimate compliance. Either size selection will deliver the desired therapy.

## Sizing: Gauntlet

**JOBST® Bella™ Strong** gauntlets are designed with overlapping sizes using wrist and palm measurements, so that patients get an optimal fit without compromising therapeutic efficacy.

Size	Palm (B)	Wrist (C)
1	15 - 17 cm	14 - 18 cm
2	17 - 19 cm	14 - 18 cm
3	17 - 19 cm	18 - 21 cm
4	19 - 21 cm	14 - 18 cm
5	19 - 21 cm	18 - 21 cm
6	21 - 23 cm	16 - 19 cm
7	21 - 23 cm	19 - 22 cm
8	23 - 26 cm	19 - 22 cm

## Sizing instructions

1. Measure palm circumference.  
This leaves 1 or 2 choices.
2. Measure wrist circumference.
3. Select gauntlet.

“

I really felt comfortable wearing it. The band at top felt much better than my original sleeve. Each day I can feel an improvement.

”




# Improved for you

## Donning made easy

It requires less force to put on the **JOBST® Bella™ Strong** Ready-to-Wear Armsleeve than those of the competition.

### Competitive Testing: Relative force to don

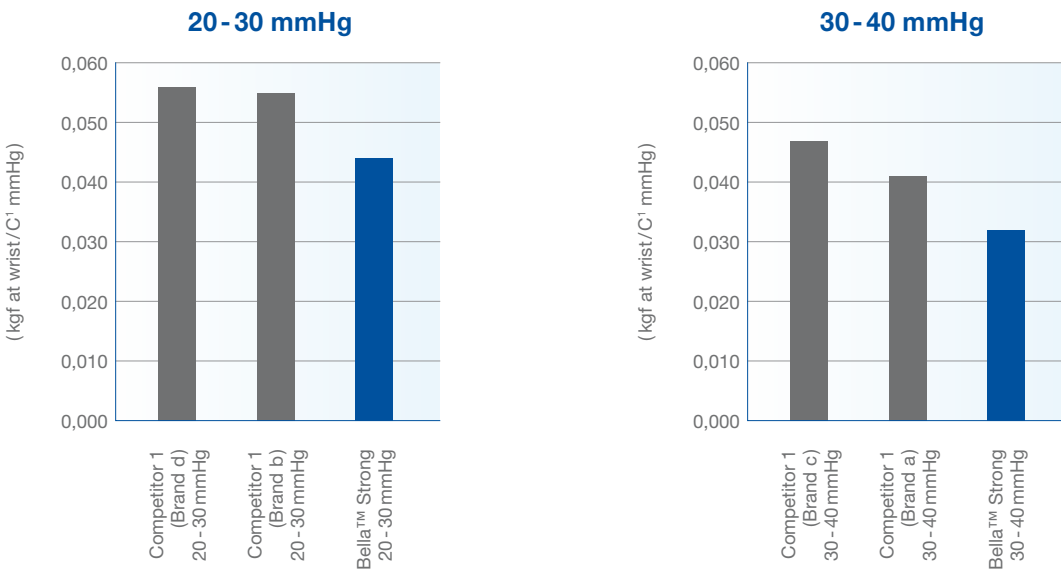

## Donning Aids

For easier donning we recommend **JOBST® DONNING GLOVES**. They create friction, provide grip and help to distribute the fabric evenly.



For easier donning we recommend using the **JOBST® DONNING AID** enclosed in package.







## Instructions for Donning

### Armsleeve


To make donning even easier each **JOBST®** Armsleeve comes with a donning aid.

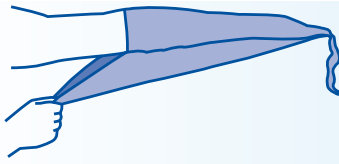


1. Pull Donning Aid onto arm.

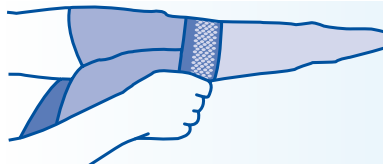


2. Fold back band on armsleeve so silicone is exposed.

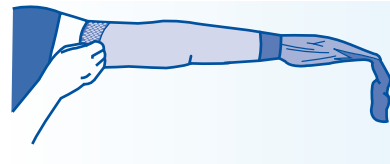

3. Pull armsleeve onto arm until wrist cuff is on or above wrist bone.

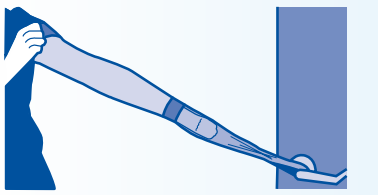


4. Place loop over door knob or other secure, stationary object. Hold the silicone band in place and slowly step away from the door – allowing the Donning Aid to be removed.

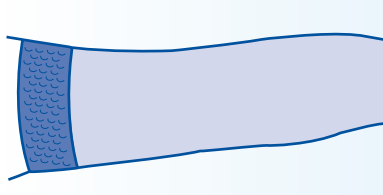

5. Flip silicone band into place using a donning glove to adjust the armsleeve as needed.

### Gauntlet

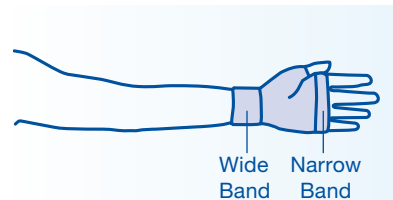

1. Slip fingers into Gauntlet starting at the end that has the widest band.
2. Position Gauntlet so that it sits comfortably around the palm and thumb.
3. The Gauntlet should be positioned so that it ends around the base of fingers.



“ No compression sleeve is going to breathe better. And I was really amazed at how well it fits after going through the washer and dryer. ”

# Improved for your world

Wide range assortment, upper extremities.

## Circular-Knit Garments

### JOBST® Ready-to-Wear Armsleeve



Compression mmHg	Beige	Black
15 - 20 mmHg	■	■
20 - 30 mmHg	■	■

### JOBST® Ready-to-Wear Gauntlet



Compression mmHg	Beige	Black
15 - 20 mmHg	■	■
20 - 30 mmHg	■	■



### JOBST® Bella™ Strong Ready-to-Wear Armsleeve



Compression mmHg	Natural	Black
15 - 20 mmHg	■	■
20 - 30 mmHg	■	■
30 - 40 mmHg	■	■

### JOBST® Bella™ Strong Ready-to-Wear Gauntlet



Compression mmHg	Natural	Black
15 - 20 mmHg	■	■
20 - 30 mmHg	■	■



# JOBST® Product Overview

## range assortment upper extremities

### Flat-Knit Garments.

#### JOBST® Elvarex® Custom Armsleeve



CCL	Compression mmHg	Beige	Black	Caramel
1	18.2 - 21.5 mmHg	■	■	■
2	25.1 - 31.5 mmHg	■	■	■
2 Forte	25.1 - 31.5 mmHg	■	■	

#### JOBST® Elvarex® Custom Glove



CCL	Compression mmHg	Beige	Black	Caramel
1	18.2 - 21.5 mmHg	■	■	■
2	25.1 - 31.5 mmHg	■	■	■
2 Forte	25.1 - 31.5 mmHg	■	■	

#### JOBST® Elvarex® Soft Custom Armsleeve



CCL	Compression mmHg	Beige	Black
1	18.2 - 21.5 mmHg	■	■
2	25.1 - 31.5 mmHg	■	■

#### JOBST® Elvarex® Soft Seamless Custom Glove / Gauntlet



CCL	Compression mmHg	Beige
1	18.2 - 21.5 mmHg	■
2	25.1 - 31.5 mmHg	■



# Bella™ Strong product overview

## JOBST® Bella™ Strong – key features

- Therapeutically effective for moderate to severe edema
- Covers unrivalled range of patient dimensions
- Armsleeve size coverage: > 80 %
- Gauntlet size coverage: > 90 %
- Ready-to-Wear
- Easy to don, armsleeve comes with donning device
- Available in armsleeve and gauntlet
- Circular-knit
- Knitted top band or 5 cm/2 inch silicone dot band
- 2 colors: natural/black

## JOBST® Bella™ Strong Gauntlet

- 15-20 and 20-30 mmHg
- 8 sizes/2 compression classes
- Comfort Flat wrist seam reduces rolling
- Extended wrist cuff pattern
- Flattened seam



Extended wrist cuff avoids pinching

Flattened seams avoid bunching

### Indications

- Edema
  - Mild to moderate Lymphedema
  - Hypertrophic scars
- Compression armsleeves can also be used for
- Phlebitis
  - Post-thrombotic syndrome
  - Vascular malformations


Silicone dot band provides secure hold, preventing curling edges and constriction



## JOBST® Bella™ Strong Armsleeve

- 15-20/20-30/30-40 mmHg
- 10 sizes per compression class
- 2 lengths: regular/long
- With or without silicone band
- Compression profile is highest at the wrist and gradient up the arm
- JOBST® Bella™ Strong sizing begins at forearm allowing for multiple size options in wrist and axilla resulting in greater size coverage
- Cotton-like feel
- Comfort Flat wristband

Supple, soft cotton-like material

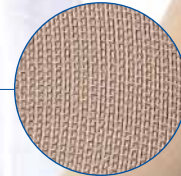

Comfort Flat wristband reduces rolling



### Contraindications for both armsleeve and gauntlet

- Arterial insufficiency
- Cutaneous infection
- Dermatitis in the acute phase
- Wet dermatosis
- Conditions where increased venous and/or lymphatic return is not desired

Please refer to the product label and/or package insert for full instructions on the safe use of these products.