HUMAN CARE

Your Life. Your Way.

User manual / Bruksanvisning /
Notice d'utilisation / Benutzerhandbuch /
Handleiding

Slings/ Selar/ Harnais / Gurte / Tilbanden
Silhouette Sling Net 25015, Silhouette Sling Trevira® 25020

Content

ENGLISH

User Manual	3
Safety Instructions	4
General Description / Intended Use	5
Technical Specifications	5
Product Label	6
Compliances and Standards	7
Suspension Alternatives	7
Hanger Bars for Human Care Lifts	8
Compatibility	9
Sling Application	10
Troubleshooting	12
Care Instructions	13

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User Manual

Important!

You must read the user manual for your device prior to usage. Keep this booklet and information handy for future use.

WARNING symbol indicates a potentially hazardous situation which, if not avoided, could result in serious injury to the user or damage to property and/or the device itself.

INFORMATION symbol indicates recommendations and information for proper, trouble-free usage.

WEIGHT CAPACITY symbol indicates maximum user weight for the product. Do not exceed this weight in any circumstance.

CE MARKING symbol indicates product meets requirements of EC directives.

RECYCLE packaging according to local regulations.

READ INSTRUCTIONS and ensure you fully understand them before using this product.

***Product modification is not permitted.
Before use, always check that the sling is
not damaged or worn.***

Human Care's products are continuously being developed and updated to provide our customers with the highest quality. We reserve the right to make design changes without prior notice.

Always make sure that you have the most recent version of the manual which is available for downloading from our website at www.humancaregroup.com.


Contact your local distributor if you have any questions about the product and its use.

Safety Instructions


The process of lifting and moving a person always involves risks. Therefore, you must read the user manuals for the lift and the lifting equipment carefully. Always make sure that the lifting equipment you are using is intended for use with the lift. As a care provider, you are responsible for the client's safety and you must know whether the client will be able to cope with the lifting process. The sling must be fitted to meet the client's specific needs. If you are in any doubt, please contact the supplier.

- Carry out a careful visual inspection of the sling every time you use it.
- Check that the right sling is being used for the client. The model, size and material must meet the client's needs.

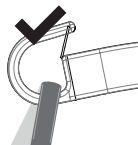Correct

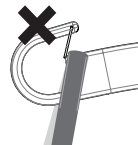Incorrect

- Check that the client is sitting securely in the sling before starting to lift him/her.
 - Check that the lift is balanced when moving, raising and/or lowering the client. If you are using a Human Care lift with telescopic hanger bar, it is important that the bar hooks are in the same position on both sides.
- Plan the moving and lifting process so that it is as safe as possible.
 - Before raising the client, check that the sling is correctly attached to the lift/hanger bar once the lift straps are stretched.
 - It is important that the hooks of the hanger bar are in the correct position, with the opening of the hook facing upwards on both sides to securely hold the straps in place. Recheck the position of the hooks once the lift straps are stretched.
 - It is important to always check that the loops are at the bottom of all hooks (see pictures) and securely placed below the locking pin on all hooks.
- Our slings have been designed and tested for indoor use with one client at a time and must only be used for transferring a person. This product is not a swing.
 - Make sure no sharp objects are in contact with the sling.

General Description / Intended Use

Silhouette slings

The Human Care Silhouette Slings are models with neck support, which are specially designed for moving people to and from narrow seating such as wheelchairs and moulded seats. Depending on the activity and the ability of the client, the arms can be inside or outside. The slings have no handles or reinforcements that may cause pressure sores. This means that it is possible to leave the sling in the chair behind the client for short periods of time.

The net slings are made of a low-friction poly-

ester net. The advantage of this material is that the sling is easy to fit and dries quickly after washing or after the client has taken a bath.


The Trevira® slings are made of a high-friction, soft cotton-like polyester. The advantage of this material is that there is less risk of the client slipping out of the sling. As a result of the friction, the sling moves relatively little in relation to the client and to the chair.

The Silhouette slings are designed for two-, three- and four-point suspension.

Technical Specifications

It is crucial to use the right size of sling for the client. To help differentiate our sizes, they are colour coded for quick identification. The colour code is found next to the sling's label in the shape of a loop. This loop can also be used to properly store or hang the sling.

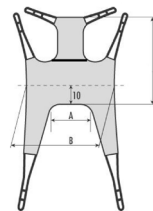Silhouette Sling

Item	Part No.:	Size	Colour code	A (cm/ inch)	B (cm/ inch)	C (cm/ inch)	Weight Capacity* (kg/lbs)	Weight Capacity (kg/lbs)
Silhouette Sling Net/ Silhouette Sling Trevira®	25015-JS/ 25020-JS	Junior small	Purple	18/7,1	62/24,4	64/25,2	15-20/33-44	250/551
Silhouette Sling Net/ Silhouette Sling Trevira®	25015-JM/ 25020-JM	Junior medium	Grey	22/8,7	67/26,4	71/28	20-35/44-77	250/551
Silhouette Sling Net/ Silhouette Sling Trevira®	25015-S/ 25020-S	Small	Red	26/10,2	77/30,3	79/31,1	30-60/66-132	250/551
Silhouette Sling Net/ Silhouette Sling Trevira®	25015-M/ 25020-M	Medium	Yellow	29/11,4	87/34,3	90/35,4	55-90/121-198	250/551
Silhouette Sling Net/ Silhouette Sling Trevira®	25015-L/ 25020-L	Large	Green	35/13,8	93/36,6	96/37,8	85-120/187-265	250/551
Silhouette Sling Net/ Silhouette Sling Trevira®	25015-XL/ 25020-XL	Extra large	Blue	40/15,7	97/38,2	106/41,7	110-140/243-309	250/551

**The weights given are only guidelines and are not a substitute for individual tests. Discrepancies can occur, tolerances ± 1 cm, 0,3".*

Slings


Product Label

The label on the sling contains the following information:

1. Human Care company logo
2. This way up arrow and marking the centre line of the sling
3. CE and MD marking
4. Product name and size
5. Picture showing the sling in use
6. Part number
7. Safe Working Load
8. "This way out, back" – outside of the sling, not in touch with the client
9. Manufacturer
10. QR-code to our webpage
11. Latex free
12. Read the manual before use
13. Visual inspection
14. Date and country of production
15. Batch number
16. Serial number
17. Washing instructions
 - a. Machine Wash Temperature
 - b. Do Not Bleach
 - c. Tumble Dry, Low Heat
 - d. Do Not Dry Clean
 - e. Never Iron
 - f. Do Not Use Fabric Softener
18. Owner name. (For example, write the client's name or another unique identifier)
19. Indicate annual inspection
20. Unique Product Code
21. Barcode/EAN

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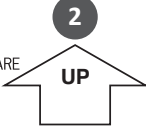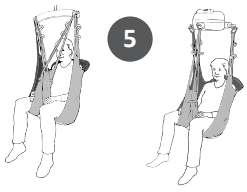
75°C / 167°F

#

Due date next periodic inspection (every 6 months), mark:

2020	1	2	3	4	5	6	7	8	9	10	11	12
2021	1	2	3	4	5	6	7	8	9	10	11	12
2022	1	2	3	4	5	6	7	8	9	10	11	12
2023	1	2	3	4	5	6	7	8	9	10	11	12

UPC

8 1608 80107 3

EAN

831599 440034

6

www.humancaregroup.com

Compliances and Standards

Human Care is an ISO 13485:2012 certified company.

This sling has been tested in accordance with ISO 10535.

The sling is CE marked in accordance with the Medical Devices Directive for Class 1 products (MDD 93/42/EEL). The CE mark is on the sling itself.

Suspension Alternatives

Our slings are designed with the premise that all people are different with varying lifting requirements, needs or desires.

To make all our clients feel as comfortable as possible in our slings, most of them can be used with 2 or 3-point suspensions. Some slings can also be used with a 4-point suspension. Below the options are displayed.

2-point suspension

Suitable for most lifting situations providing a good seating position. 2-point suspension is available with floor lifts, fixed and portable ceiling lifts.



Two-point suspension with crossed leg support:
Thread the strap on one leg section through the strap on the other leg section before attaching them to the hanger bar.

Two-point suspension with overlapping leg support:
This method is suitable for lifting clients with amputated legs, for example. Extension straps can be attached to the lifting straps to prevent the client from leaning too far back in the seated position.

Two-point suspension with separate leg support:
Put the client's legs in each leg section and attach the lifting straps to the hanger bar. Note: This method increases the risk of the client sliding out of the sling

An individual assessment of the client's needs must always be carried out before choosing an alternative connection method.

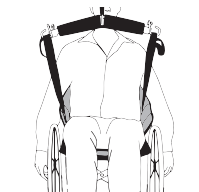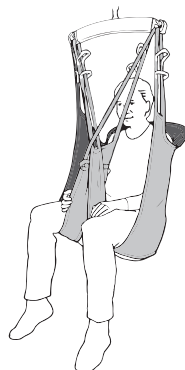
Slings

3-point suspension


Offers an elongated sitting position and more even weight distribution. The 3-point suspension is suitable especially for clients carrying an increased risk of pressure sores or who need additional space for the upper body.

Thread the left leg section through the loop on the right leg section. Attach the leg section lifting straps to the leg hooks and the back section lifting straps to the lift.

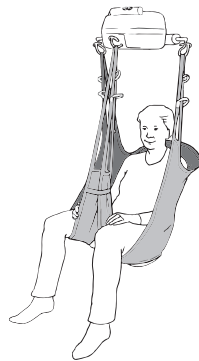
4-point suspension


A spacious sitting position making it suitable for larger clients or when using a stretcher. 4-point suspension is available with both fixed and portable ceiling lifts or floor lifts.

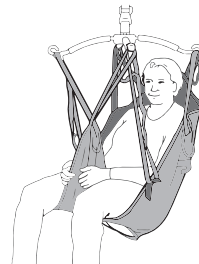
Hanger Bars for Human Care Lifts

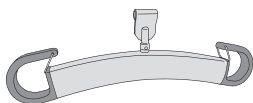HeliQ
(7275x)

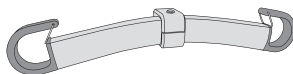Roomer S / Altair
(55751 45 cm, 17,7" / 55764 55 cm, 21,7")

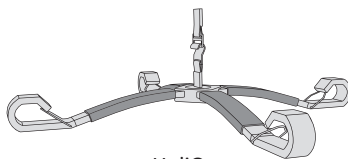HeliQ
(70800)

Accessory	Description	Part no.
2-point hanger bar	For stationary lift HeliQ. 3 widths:	35 cm/13,8": 72754 45 cm/17,7": 72755 60 cm/23,6": 72757
2-point hanger bar	For portable lifts Roomer S / Altair, 2 widths:	45 cm/17,7": 55751 55 cm/21,7": 55764
4-point hanger bar	For stationary lift HeliQ	70800
4-point suspension	4-point suspension accessory (49 cm). Enables 4-point lifting with the Altair and Roomer S lifts. The suspension bar is placed on the lift's single leg hook.	55619

Compatibility

Human Care Group slings are designed for use with all Human Care lifts and hanger bars. The products are also compatible with products of other suppliers utilizing the same attachment method for securing slings to a hanger bar and lift. However, due to the wide range of products worldwide, Human Care cannot be responsible for any errors that may result from improper use or application of combinations of slings, hanger bars and lifts of other suppliers. The use of the combined system is at the prescriber's own risk and liability.

To ensure the safety of the client and caregiver, Human Care recommends the following minimum requirements:

1. When using slings with lifts and hanger bars, made by manufacturers other than Human Care, the combined system needs to be designed with the same sling attachment method for the hanger bar of floor lifts and ceiling lifts, meaning loop style slings are combined with a loop style hanger bar, clip style slings are combined with a clip style hanger bar.
2. The involved patient lifts and sling must be CE marked in accordance with the medical device directive, MDD 93/42/EEG / upcoming MDR 2017/75 and comply

with standard requirements for patient lifts as stated in EN ISO 10535.


3. The caregiver must follow the manufacturer's instruction manuals and recommendations with respect to the use, care, client, caregiver, cleaning, and inspection of the slings and lifts.
4. Different products on the assembled lift system: i.e., lift, hanger bar, sling, scale and other accessories may have different maximum loads. It is always the lowest maximum allowable load indicated for the respective products and the system that applies for the entire system. Always check the markings on all individual products in the system.
5. Individual risk assessment, including a practical test of the system and validation with the specific combination used are required to affirm that the size of the sling is correct for the intended use and compatible to the width and design of the hanger bar. The combined system needs to be documented by the competent assessor.

If there are any questions, please contact your local Human Care office or distributor.


Sling Application

Lifting the client from a lying position

1. Lay the client on his/her side. Position the sling so that the centre of the sling lines up with the client's spine. Take note of the marking on the sling "This way out, back". The bottom part of the sling must be directly under the client's coccyx. Roll the client over onto his/her other side and pull the sling through. Lay the client on his/her back. If the bed is height-adjustable, it is a good idea to raise the head of the bed so that the client is in a semi-sitting position.

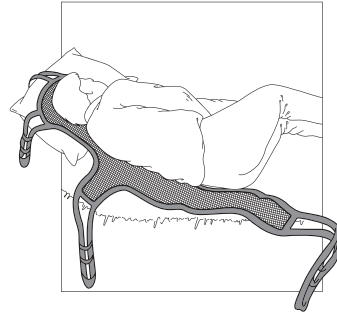
2. Pull the leg sections through between the client's legs. Lower the lift to the height of the client's stomach. When a three-point suspension is used, you should thread the left leg section through the loop on the outside of the right leg section. Attach the leg section lifting straps to the leg hook on the lift. Attach the back section lifting straps to the lift. Hold the lift away from the client's face with one hand as you raise it.

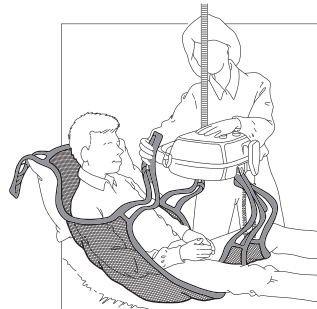
Alternative: The sling can also be used for two- or four-point suspension. In this case the loop on the right leg section is not used and the leg sections are crossed.

3. Attach the neck support lifting straps to the lift bar. Pause in this position and check that all the lifting straps are securely attached to the hooks.

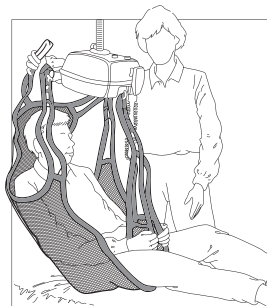- Continue lifting the client. It may be necessary to lower the bed slightly. Raise the client's legs slightly while turning him/her away from the bed. It is a good idea to put your arm around the client while you are moving him/her to create a greater sense of security.

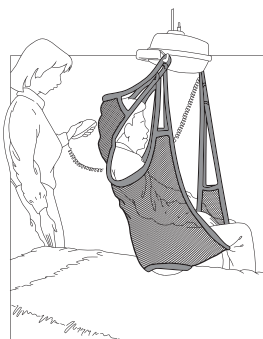
Lifting the client from a sitting position

- Pull the sling down behind the client's back. Take note of the marking on the sling (outside up). The bottom edge of the sling must be behind the client's coccyx. Pull the leg sections through to the front.


- Pull the leg sections up between the client's legs. Lower the lift to the height of the clients's stomach. When a three-point suspension is used, you should thread the left leg section through the loop on the outside of the right leg section. Attach the leg section lifting straps to the leg hook on the lift. Turn the lift around and attach the back section lifting straps onto the lift bar. Hold the lift away from the client's face with one hand as you raise the lift. Alternative: The sling can also be used for two- or four-point suspension. In this case the strap on the right leg section is not used. Instead the leg sections are crossed.

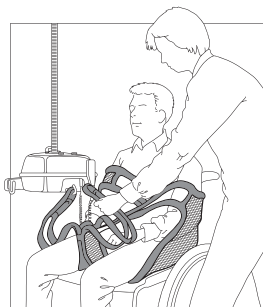
Slings

3. Attach the neck support lifting straps to the hanger bar. Pause in this position and check that all the lifting straps are securely attached to the hooks.

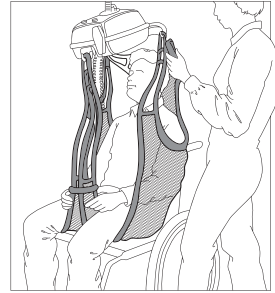4. Continue lifting the client. It is a good idea to put your arm around the client while you are moving him/her to create a greater sense of security.

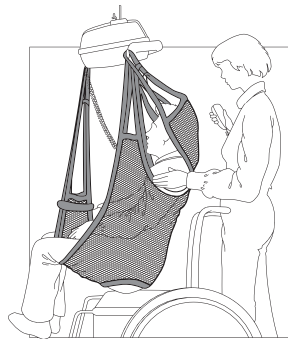Troubleshooting

Description	Action
The client is at risk of sliding out	<ul style="list-style-type: none">• The sling may be too large.• The sling is not correctly fitted- the length of the loops should be stretched at the same time.
The client is unbalanced when lifted	<ul style="list-style-type: none">• The instructions on the individual test log have not been followed or a new assessment needs to be carried out.
The sling is too small	<ul style="list-style-type: none">• This can be uncomfortable for the client. The sling can cut into the client's flesh and/or his/her back may not have enough support.

Note: It is important that the right sling is used for each client and that it is individually tested to ensure that it meets the client's need and the transfer situation.

Care Instructions

Check for wear and tear

Check the sling for wear and tear. It is important to check the fabric, the edging, the seams and the lifting straps.

Damaged lifting equipment must not be used.

Cleaning

Human Care slings are intended for individual use. To avoid infections or transmission of bacteria, it is important to clean the slings properly. All slings should be machine washed and dried as soon as possible. Slings should always be cleaned after use in a wet situation.

Always refer to the label of the sling before washing.

Periodic inspection

A periodic inspection must be carried out at least every 6 months. See Human Care's peri-

odic inspection log which can be found on our website at www.humancaregroup.com

Expected sling service life

The sling is expected to function for approximately 150-180 laundering cycles, depending on the intensity of use, maximum load applied during use, and adherence to correct laundering procedures. This equates to approximately 3 years of laundering once per week or 5 years of bi-weekly laundering. It is recommended that the sling is washed separately using a hospital grade laundry bag. If the product label is no longer legible, the sling should be discarded.

Recycling

Once the product has reached the end of its life, it must be disposed as residual waste.

Overview of Sling Material Care		Polyester	Polyester Net	Trevira	Polyester/Polyester Net	Polyester/Cotton	Ventilated Polyester	Cotton	Polyester/Triac	Coated Nylon
	Machine Wash Maximum Temperature 60°C (140°F) (slings with padding)	●	●							●
	Machine Wash Maximum Temperature 75°C (167°F)	●	●	●	●		●		●	
	Tumble Dry Low, Maximum Temperature 30°C - 54°C (86°F - 130°F)									●
	Tumble Dry Low, Maximum Temperature 60°C (140°F)	●	●	●	●		●		●	
	Do Not Tumble Dry					●		●		
	Do Not Iron	●	●	●	●	●	●	●	●	●
	Do Not Dry Clean	●	●	●	●	●	●	●	●	●
	Do Not Bleach	●	●	●	●	●	●	●	●	●
	Do Not Use Fabric Softener	●	●	●	●	●	●	●	●	●
	Do Not Wash					●		●		
	Dispose of After Use					●		●		
	Disinfection Permitted (see list of approved disinfectants)	●	●		●		●		●	

Slings

Disinfection Procedure

Human Care slings are intended to be used with one patient at a time. To avoid possible transmission of infectious diseases, it is important to clean the sling if it has become soiled or if it is to be used by a different patient.

Approved Disinfectants	Solutions to Avoid
Chemical Disinfection: 70% Ethanol, 45% isopropanol or similar	 Steam sterilization
Thermal Disinfection: Wash 65°C (149°F) for a minimum of 10 minutes ¹	 Iodine based disinfectants
Thermal Disinfection: Wash 71°C (160°F) for a minimum of 3 minutes ¹	 Gas sterilization
Note: Add 8 minutes to both thermal disinfection options to ensure heat has penetrated the fabric. Finish with a rinse cycle. ¹	

¹ Ref: HSG (95) 18: Hospital Laundry Arrangements for Used and Infected Linen, Department of Health, NHS (National Health Service, England)

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