

ALENTI

MULTIPURPOSE
HYGIENE CHAIR

ARJO



WHAT THE ALENTI DOES

With the *Aleuti* lift hygiene chair, one nurse can safely bring residents who are able to sit up, from bed to ARJO hygiene equipment and back again, without manually lifting, minimizing stress and strain on the caregiver. A moveable back-rest and safety belt allows for transfers from either side of the bed. The *Aleuti* is battery powered for height adjustment via a hand-set. A removable and rechargeable battery makes the lift hygiene chair ready for use at all times. The *Aleuti* operates as an integrated system with all lengths of ARJO height adjustable recumbent bathing systems (*Rhapsody* and *Primo*) and with the ARJO *Prelude* shower cabinet as well as the ARJO *Century Bath*. The *Aleuti* can also be used as a shower chair, allowing the nurse to work

in an ergonomically correct posture, to prevent static load injury. An optional electronic scale facilitates the weighing of residents.

The *Aleuti* is a versatile mobile lifting aid which can transfer residents directly from their bed or wheelchair. Use of the *Aleuti* not only eliminates manual lifting, but also makes it possible to reduce the number of transfers. No intermediate transfer movements are necessary from the moment the resident is safely seated on the *Aleuti* at collection for bathing or showering, to the moment he or she is transferred back to bed.

The *Aleuti* is delivered fully assembled and ready for immediate use.



DESIGN

Leading by innovation and design has been a characteristic of ARJO over the years. ARJO *Aleuti* has received a number of different design awards. Among others, the Industrie Forum Design award and Excellent Swedish Design.



ERGO GRIPS

Ergonomically designed handlebar grips allow for easy movement of the *Aleuti* at the bedside, the bath or around corners.



FEEL BACK SEAT DESIGN

Those residents with a greater degree of mobility are able to seat themselves without difficulty due to the ergonomically designed seat and chassis design, which allow the resident to stand close to the lift, which make sitting down and standing up much easier.





KEEP GOING WITH 2XPOWER

Alenti is always delivered with two high capacity batteries. This ensures the *Alenti* is always ready for use.



ON-BOARD SCALE

The *Alenti* can also be used to weigh residents. It utilizes a highly accurate on-board electronic scale. Weighing, a difficult task, can now be easily performed on the way to the bath or shower, without the residents clothing, for more accurate weights.



EASY CLEAN

Its one-piece construction and smooth surface, made out of smooth closed cell polyurethane rubber, makes it hygienic and easy to clean and disinfect.



HANDY BRAKE

Raising, lowering and central braking functions are controlled with a hand control, enabling the caregiver to stay close to the resident.



SAFETY BELT

The safety belt ensures that the resident remains in a safe and secure position while on the lift and prevents the resident from sliding on the lift while bathing.




FREE WHEELING

Specially designed wheels makes the *Alenti* more maneuverable and makes turning easier.



PRODUCT SPECIFICATIONS

Space requirement

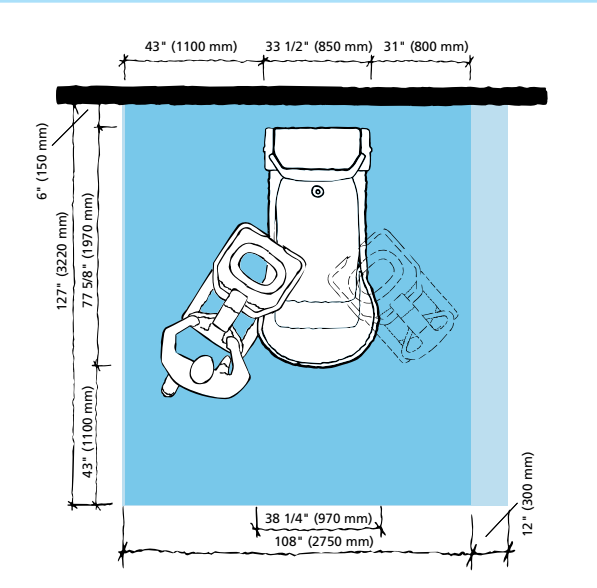

Resident Gallery

Suitable for Albert, Barbara and Carl



- A Albert** who is ambulatory and independent
- B Barbara**, who can support herself to some degree
- C Carl**, who sits in a wheelchair and has little capacity to support himself

Please contact ARJO for further information on the Resident Gallery™.

Product information

Safe working load	352.7 lbs. (160 kg)
Lifting stroke range	23" (585 mm)
Width	23.81" (605 mm)
Max length, highest pos.	35.43" (900 mm)
Supply voltage	24 V DC

Visual brake indicator

Operated both from hand set and/or pillar mounted controls

Handset (protected to class IPX6)

Safety belt

Low-friction castors

Seat is made of a special blend of polyurethane resin for warmth, comfort, and superior infection control

Alenti integrates with the following ARJO products:

Century

Rhapsody / Primo

Prelude Shower Cabinet

SARA 3000

SARA Plus

Alenti accessories:

Integrated scale

Bedpan holder

Bedpan with lid

Leg rest for pedicure

Approvals

UL, cUL, TÜV

Use only ARJO Designed Parts, i.e. parts specifically designed for the purpose, on equipment and products supplied by ARJO. Our Conditions of Sale specify that ARJO is not liable for injuries attributable to the use of inadequate parts. Continuous development is part of our policy, and we therefore reserve the right to change specifications without notice.

Alenti, Century Bath, Encore, Prelude, Primo, Rhapsody and Resident Gallery are trademarks belonging to ARJO International AG

ARJO

...with people in mind.

ARJO, Inc.
 50 N. Gary Avenue
 Roselle, Illinois 60172
 Phone: 800-323-1245
 Fax: 888-594-2756
www.arjo.com

MEMBER OF THE GETINGE GROUP