

Altus Pro System

True Low Air Loss Mattress System

The Altus Pro System is a true low air loss mattress system suitable for medium and high risk pressure ulcer treatment and 24 hour pressure area care. It has been designed to reduce interface pressure and aid in the treatment and prevention of pressure ulcers/sores; specifically patients with visible pressure ulcers/sores falling under Stages II – IV.


The mattress and blower system automatically maintain four different pressure zones to reduce the incidence of pressure ulcers while optimizing patient comfort (softer head and heel sections, and firmer torso and leg sections). Ideal for improving circulation, lowering body temperature and promoting a dryer environment.

Great for nursing home, long term care and home care settings as prescribed by a physician.




Altus Pro System

FEATURES

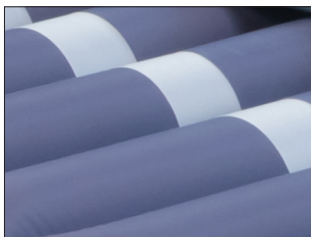

Altus Pro Mattress

- 20 rows of 10" nylon/PU cells
- Provides a continuous pressure relief environment while helping to reduce body temperature and moisture by allowing air to diffuse and circulate through the laser holes of the mattress surface.
- Mattress base has an integrated loop bridge to hold the air cells in place.
- Snap cell design allows for easy replacement of individual cell.

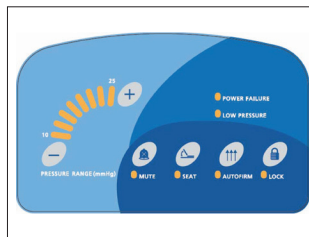

Altus Pro Pump

- Lightweight compact design attaches easily to bed ends.
- True Low Air Loss Blower produces over 1000 LPM in air flow.
- Fully inflates mattress in under two minutes.
- 8 levels of pressure settings allowing pressure to be customized to patient's weight.
- Auto Firm Function quickly provides uniform firmness.
- Seat Function increases pressure when patient is in the sitting position to prevent bottoming-out and for increased patient comfort.
- Visible and audible alarms with Alarm Mute Function.
- Emergency CPR Deflation through quick connector within 15 seconds.



Altus Pro Cover

- 4 way stretch polyester PU transfer cover with waterfall zipper for added protection.
- Non-skid surface on the base of cover to keep mattress secure when the bed is raised at either end.
- Fluid resistant
- Vapor permeable
- Anti-microbial

SPECIFICATIONS

Pump

| | |
|--------------------|--------------------------------------|
| Power: | 120VAC/60Hz/250W |
| Input Frequency: | 60 Hz |
| Air Volume: | Greater than 1,000 liters per minute |
| Pressure Range: | 10mmHg ~ 25mmHg |
| Weight: | 10.5 lbs. |
| Size: | 10.7" x 6.3" x 12.2" |
| Power Cord Length: | 14.8 ft. |
| Warranty: | 2-year |
| UPC: | 092237602966 |

Mattress

| | |
|----------------------|------------------------------------|
| Size: | 80" x 35" x 10" |
| Number of Cell Rows: | 20 rows of 10" high nylon/PU cells |
| Mattress Weight: | 22 lbs. |
| Cover Material: | Polyester PU Transfer |
| Cell Material: | Nylon 70D/PU |
| Base Material: | Nylon 840/PU |
| Max Weight Capacity: | 425 lbs. |
| Warranty: | 2-year |

Billing

| | |
|-------------|-------|
| HCPSC Code: | E0277 |
|-------------|-------|

⚠ WARNING

Whenever bed rails are used as part of any hospital or home care bed system, the potential for bed rail entrapment exists. To reduce the risk, a patient must always be assessed regarding their risk of entrapment before any product is used. Routine patient monitoring is also required once a product is in use. To learn more about the risk of bed rail entrapment, contact your Roscoe Medical sales representative for a copy of our "Bed Rail Safety Guide" or visit the U.S. FDA's website at <http://www.fda.gov> and search for "Bed Rail Entrapment".

| Item # | Description | Qty |
|----------|---|------|
| APM-TLAL | Altus Pro True Low Air Loss Mattress System | each |

* Providers are responsible for determining appropriate billing codes for Medicare Program claims. This information is included for your convenience and is not intended, nor should it be considered billing advice.

