



# Domus 1

## Alternating Pressure Redistribution System

- Suitable for patients with low risk of pressure ulcers
- Medical grade PVC pad for better hygiene and safety
- An ideal transition mattress for homecare



Care for a Healthy Life


**APEX**

# Domus 1

## Alternating Pressure Redistribution System



**Anti-Sliding Extension Flap**  
Can be tucked under the standard mattress to prevent pad from sliding

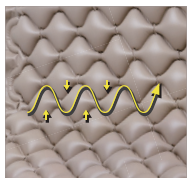


**Medical Grade PVC Pad**  
Effortless for cleaning and fire retardant for better safety



**Adjustable Hangers**  
Can secure control unit to fit most of beds



**Weight Indicator Reference**  
To facilitate optimized pressure setting

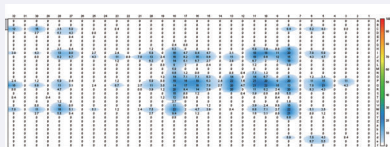

**Pressure Redistribution During Alternating**

Pressure mapping reveals that bladder A & B take turns to inflate and deflate can significantly relieve the body interface pressure

### Specification

Pump	Dimension	9.8" x 4.9" x 3.3" / 25 x 12.5 x 8.5 cm
	Weight	2.4 lb / 1.1 Kg
	Case Material	Flame retardant ABS
	Supply Voltage	AC100-120V 60Hz, 0.1A
	Operating Cycle	10 min
Mattress	Dimension	2.5" Bubble Pad Overlay
	Dimension	77.2" x 35.4" x 2.5" / 196 x 90 x 6.4 cm
	Cell Height	2.5" x 130 pcs of bubbles
	Weight	4.6 lb / 2.1 Kg
	Material	Medical - grade PVC
	Maximum Patient Weight	220 lb / 100 Kg
	Flame Retardant Standards	CAL 117
HCPCS code		E0181

*Specifications are subject to change for improvement without notice*

### Ordering Information (For North America Only)

	Domus 1
System	PP01101
Mattress	PM03401
Pump	PP01001