

MATTRESS LATERAL ROTATION



COMBINES 3 THERAPIES

- LATERAL ROTATION
- ALTERNATING PRESSURE
- LOW AIR LOSS
TO EFFECTIVELY TREAT PRESSURE ULCERS.



LATERAL ROTATION

Rotates the patient from side to side to provide proper pressure redistribution.



ALTERNATING PRESSURE

Redistributes pressure throughout the length of mattress.



LOW AIR LOSS

Treats bed sores by creating a dry cool environment.




Customized Lateral Rotation



Side bolsters for optimum patient safety reduces the risk of falling.

Mid-way transfer point allows for easy patient egress and ingress

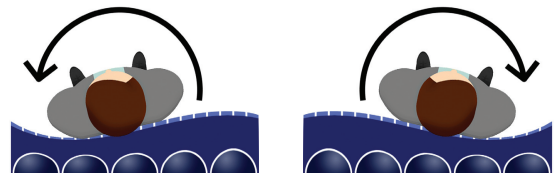

MATTRESS LATERAL ROTATION



PUMP



Alternating Pressure and Static Function



Comfort Level Adjustable



Low Air Loss with on/off function



5 Level Cycle Time Adjustable



- 5 Level Rotation Time Adjustable
- Rotation On/Off Function
- Panel Lock Function
- Power Failure Alert
- Low Pressure Alert
- Audible Alarm Mute Function
- Quick-change Air Filter

PRESSURE RANGE

Digital readout adjusts pressure range from 30 - 60mmhg

ALTERNATING CYCLE TIME

5, 8, 10, 15, 20 minutes selectable

ROTATION CYCLE TIME

15, 30, 45, 60, 90 minutes selectable

AIR OUTPUT

10 LPM

PUMP DIMENSIONS

11" x 12" x 5.9"

MATTRESS



Alternating Pressure



Low Air Loss



Lateral Rotation



Pillow Function - Top 3 cells remain static for ultimate patient comfort



Adjustable Straps secures mattress to bed frame



LOW AIR LOSS

8 Ventilated Cells

CELL NUMBER

20 Easy Replaceable Air Cells

CELL MATERIAL

TPU

COVER MATERIAL

2 Way Stretch PU Fluid Resistant Cover

BASE MATERIAL

Woven with coated non-slip material

CPR

Pull strap

DIMENSIONS

80" L x 36" W x 8" H with 2" Side Bolsters

WEIGHT CAPACITY

440 Lbs.



Safe yet effective Angle of Rotation: 30 degrees



MODEL # **LRM36**