marcy.

There for you from the start.















I am your anterior walking aid, and will assist you in your very first attempts


to walk. I can hardly wait to start off with you! I will support your balance and simultaneously train your torso musculature. You can adapt my wheelbase to suit: if you adjust it so that it is somewhat shorter, I become highly manoeuvrable and you can move me better. I give you more security and support with a long wheelbase. If you find it hard to walk in a straight line, I am happy to assist you with my direction locks.



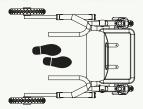



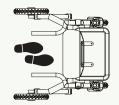

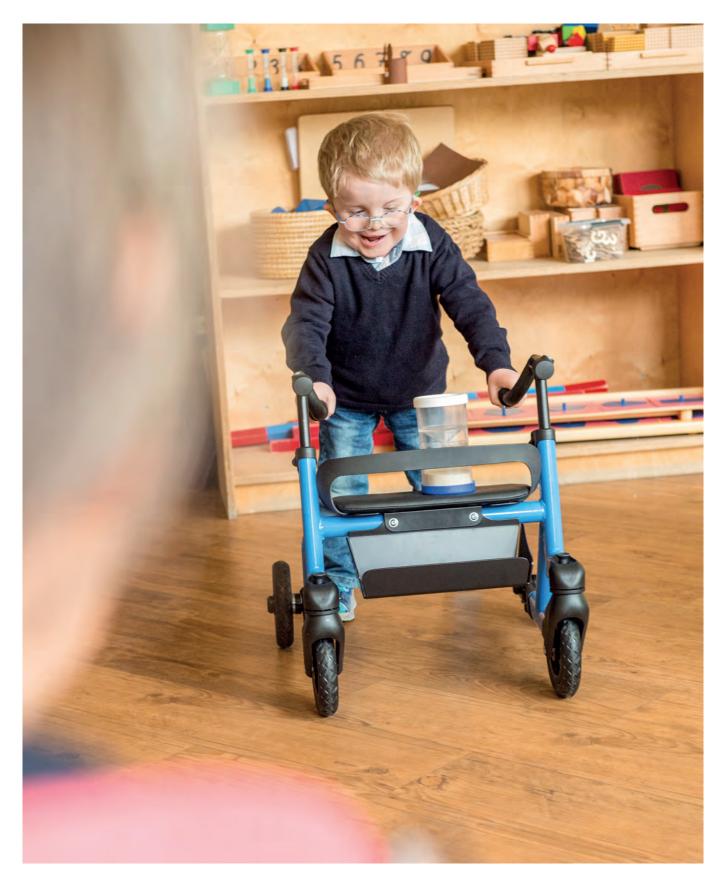



marcy.





Optional parking handbrake provide you with more safety and control - also in combination with the forearm supports.



Frame colours Size 0



Frame colours Size 1













| Data | Size 0 | Size 1 | Size 2 | Size 3 |
|-------------------------------------|---------------|-------------|--------------|-------------|
| Handle height | 11.7"-18.3" | 14.8"-23" | 20.6"-29.25" | 24.5"-34.3" |
| Height of forearm rests | 15.9"-22.6" | 19.1"-27.3" | 24.5"-33.1" | 28.4"-38.2" |
| Handle spacing | 7"-12.4" | 8.5"-14" | 11.7"-17.1" | 14"-19.5" |
| Seat height | 9.75" | 11.7" | 17.5" | 17.5" |
| Total width | 17.1" | 18.7" | 24.1" | 27.6" |
| Total length | 14.2" | 16.5" | 22.4" | 23.7" |
| Total length (telescopic wheelbase) | 15.2"-18.3" | 17.5"-21.4" | 23.4"-28.8" | 25.3"-30.8" |
| max. load | 44 lb. | 66.1 lb. | 110 lb. | 143.3 lb. |
| Weight | 9.2 lb. | 11 lb. | 16.3 lb. | 20 lb. |
| HMV No. | 10.46.02.3049 | | | |