



slings

molift[®]
designed for life



Our Philosophy

The sling is like a second skin for the user. It does not matter how good the lifter is, if the lifting sling is not correct. That is why we focus on developing the most flexible and comfortable slings for both users and care-givers. Our goal is to provide comfort, safety and dignity for the user during the lifting and transfer. For this reason we have developed a 4-point suspension as a standard on our lifters. This makes the lift very stable, and the user will not swing from side to side.

The 4-point suspension will also allow a free and open sitting position, especially around the head area without any risk of bumping the head into the suspension. The user will not feel squeezed as they would in a 2-point suspension, but comfortable and safe. With the combination of sling and suspension, the user will slide into a natural sitting position.

The slings are designed to distribute weight and pressure evenly.

Molift offers a wide, flexible selection of slings for all kind of transfers. The slings vary in models and sizes. Because Molift focuses on comfort and safety we strongly recommend an individual assessment. A colour code indicates the sling size, which makes it easier and quicker to choose the correct sling size. We design and manufacture our own slings in Norway to the same high standards as the Molift lifters, and all our slings are approved in accordance with the NS-EN ISO 10535. Molift recommend the use of Molift slings in combination with Molift lifters.

Molift Easy Sling

Molift Easy Sling is suitable for most lifting situations. It is available with and without head support.

- Padded for more comfort

The body and leg parts have a firm padding to make it more comfortable for the user. It distributes weight and pressure evenly, and makes it less likely that the sling will indent, pinch or otherwise be disagreeable against the user's skin.

- Just the easiest sling to use

A simple, unique folding technique makes it possible to use one hand for the most critical phase of fitting the sling around the user, leaving the other hand free to support and safeguard the user.



Molift Easy Sling. Information and part numbers

	without head support	with head support
XS	Part No. 3025000	Part No. 3026000
S	Part No. 3025100	Part No. 3026100
M	Part No. 3025200	Part No. 3026200
L	Part No. 3025300	Part No. 3026300
XL	Part No. 3025400	Part No. 3026400
XXL	Part No. 3025500	Part No. 3026500



Use positioning loops on the leg parts to achieve a reclined sitting position.



Open and comfortable sitting position.



We recommend Molift Easy Sling with head support when lifting from a lying position.



Anatomic design respecting the natural curve of the back, makes the sling comfortable and prevents the body from gliding.



Molift Easy Toilet Sling

Lifting and transferring users to and from the toilet can be the most difficult and least dignified lift. Compared to a regular sling the lifting points are positioned differently - from the thighs and seat to the hollow of the knees and high on the chest to obtain the necessary opening around the seat area to adjust clothing.

- Multiple positioning loops

Adjustable leg and shoulder straps to achieve a 90 degree sitting position which ensures a correct toilet positioning.

- Twin belt system

Two buckles on the waist band secure an even distribution of the pressure and the ability to tighten correctly for increased safety and stability.

- Flexible leg parts

The hinged leg parts facilitate the adjustment of clothing.




Molift Easy Toilet Sling. Information and part numbers

	without head support	with head support
XXS	Part No. 3032050	Part No. 3033050
XS	Part No. 3032000	Part No. 3033000
S	Part No. 3032100	Part No. 3033100
M	Part No. 3032200	Part No. 3033200
L	Part No. 3032300	Part No. 3033300
XL	Part No. 3032400	Part No. 3033400
XXL	Part No. 3032500	Part No. 3033500



Molift Easy Toilet Sling with head support.

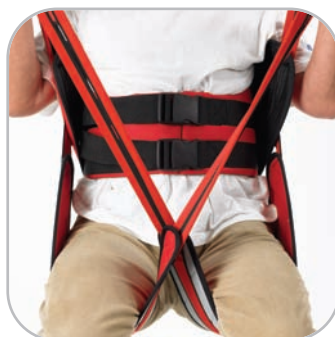

Twin belt system.



Good back support and flexible leg parts.



Multiple positioning loops.

Molift Easy Amputee Sling

Molift Easy Amputee Sling is designed for safe, comfortable transfers of users with single or bi-lateral amputation. Compared to a regular Molift Easy Sling the opening around the seat area is reduced, and also the leg parts are custom made for amputees.

- Twin belt system

Two buckles on the waist band secure an even distribution of pressure and the ability to tighten properly increasing safety and stability.

- Reduced aperture

Molift Easy Amputee Sling is suitable for bi-lateral amputees.



Molift Easy Amputee Sling. Information and part numbers

	without head support	with head support
XXS	Part No. 3027050	Part No. 3028050
XS	Part No. 3027000	Part No. 3028000
S	Part No. 3027100	Part No. 3028100
M	Part No. 3027200	Part No. 3028200
L	Part No. 3027300	Part No. 3028300
XL	Part No. 3027400	Part No. 3028400
XXL	Part No. 3027500	Part No. 3028500



Without head support.



Reduced aperture and loops to prevent leg parts from sliding apart.



Molift Pediatric Slings

Molift Pediatric Slings are designed especially for children. They are available with and without head support.

- Same wide range as adult slings

All Molift slings are available in pediatric sizes.

- Tutti-Frutti colours

Molift Easy Pediatric Sling has cheerful colors to distinguish the slings from the regular range.



Molift Easy Pediatric Sling. Information and part numbers

	without head support	with head support
XXS	Part No. 3005050	Part No. 3006050
XS	Part No. 3005000	Part No. 3006000

Molift Easy Toilet Sling. Information and part numbers

	without head support	with head support
XXS	Part No. 3032050	Part No. 3033050
XS	Part No. 3032000	Part No. 3033000

Molift Basic Bathing Sling. Information and part numbers

	without head support	with head support
XXS	Part No. 3023050	Part No. 3024050
XS	Part No. 3023000	Part No. 3024000



Molift Easy Pediatric Sling with head support and cushion.



Molift Easy Toilet Sling.



Molift Basic Bathing Sling.

Molift Basic Sling

Molift Basic Sling is the same shape and size as Molift Easy Sling, but without padding. Molift Basic Sling should be chosen when the padding makes fitting of the sling more difficult, such as lifting and transfer from moulded seats in wheelchairs or comfort chairs.

- Non-padded alternative

Molift Basic Sling is appropriate for most lifting and transferring situations.

- Maintains the positioning abilities

Same possibility to achieve an upright sitting position as with Molift Easy Sling, as well as the adjustment for shortening the leg straps to achieve a more reclined position.



Molift Basic Sling with head support.

Molift Basic Sling. Information and part numbers

without head support

XXS	Part No. 3021050
XS	Part No. 3021000
S	Part No. 3021100
M	Part No. 3021200
L	Part No. 3021300
XL	Part No. 3021400
XXL	Part No. 3021500

with head support

Part No. 3022050
Part No. 3022000
Part No. 3022100
Part No. 3022200
Part No. 3022300
Part No. 3022400
Part No. 3022500



4-point suspension ensures an open and comfortable sitting position.



Shoulder loop to lock sitting position.



Molift Basic Sling may make fitting easier in a comfort chair.

Molift Basic Bathing Sling

Molift Basic Bathing Sling is the same design and shape as the Molift Basic Sling.

- Dries quickly

Molift Basic Bathing Sling is made from polyester mesh net which is comfortable next to the skin.

- Allows both reclined and upright position

Molift Basic Bathing Sling is easily adjustable due to multiple positioning loops.



Molift Basic Bathing Sling. Information and part numbers

	without head support	with head support
XXS	Part No. 3023050	Part No. 3024050
XS	Part No. 3023000	Part No. 3024000
S	Part No. 3023100	Part No. 3024100
M	Part No. 3023200	Part No. 3024200
L	Part No. 3023300	Part No. 3024300
XL	Part No. 3023400	Part No. 3024400
XXL	Part No. 3023500	Part No. 3024500



Molift Basic Bathing Sling with head support.

Molift Hammock Sling

Molift Hammock Sling enables a comfortable, fully supported sitting position. It adapts to the contours of the body, and the user can be left in the sling after the transfer, for example, to a comfort chair or a wheelchair. This kind of sling is often custom made.

- Breathing fabric

The sling can be left under the user for a longer period. It is made of a polyester net which allows body moisture and heat to come through.

- Lining for extra comfort

Molift Hammock Sling can be delivered with two different soft, skin friendly detachable linings as an option, either in wool or in fleece.



Molift Hammock Sling. Information and part numbers

One size	Part No. 3047000 - according to specifications
One size	Part No. 3047001 - wool lining according to specifications
One size	Part No. 3047002 - fleece lining according to specifications



The sling can be left under the user.



Molift Hammock Sling with fleece lining.



Molift Patient Specific Sling

The Molift Patient Specific Sling is a disposable sling which is the same shape as our regular Molift Slings.

- The ultimate sling for infection control

The Molift Patient Specific Sling is used by the individual patient; thus helping to reduce cross infection. The contaminated sling should be disposed of in accordance with current infection control procedures. If the sling is subjected to any type of fluid, the inner white fabric will change color to alert the caregiver that the sling should be replaced immediately with a new one. It is not designed to be washed or re-used.

- Padded for more comfort

The Molift Patient Specific Sling is a comfortable, padded high quality sling with adjustable head support. It is manufactured from an environmentally friendly material.

Molift Patient Specific Slings can also be used as an alternative to laundering slings. They are sold in packs of 10.



Molift Basic Bathing Sling with head support.



Molift Patient Specific Sling. Part numbers

	with head support
S	Part No. 3050100
M	Part No. 3050200
L	Part No. 3050300
SWL	660 lbs

Molift Fabric Stretch



Molift Fabric Stretch is designed for lifting from a lying position. It is a functional stretcher especially designed for lifting patients who must be kept in a stable position. Molift Fabric Stretch must be used with a special 8-point suspension.

- The comfortable horizontal lift

Molift Fabric Stretch is a comfortable alternative to stretchers and flatlifters.



Molift Fabric Stretch. Part numbers

352 lbs Part No. 3040001

660 lbs Part No. 3040002

Molift Stretcher



The Molift Stretcher is a spinal-scoop stretcher that is perfect for lifting users with suspected fractures. It ensures safe patient handling from the floor or bed. The stretcher locks at the top and bottom, and can be split in two parts.

The Molift Stretcher can be used in X-Ray and CT scanning.

- Easy and user friendly

Molift Stretcher has a high lifting capacity of 300 kg. It can be used directly on the standard Molift 4-point suspension on Molift Partner 205, 255 and Molift Nomad. Security straps can be ordered separately. The stretcher is also available with reels that automatically retract the lifting straps when not in use.



Molift Stretcher.

Measurements: 78 x 25 x 126 in (L x W x H)

Total weight: 32 lbs, half: 16 lbs

Material: Plastic composite and stainless steel locks.

Part No. 2150103

SWL: 660 lbs

Molift Ambulating Vest

Molift Ambulating Vest is specially developed for rehabilitation. It is also suitable for toilet situations.

-Safe rehabilitation

Molift Ambulating Vest has adapted features from professional mountain climbing straps to make it extra safe and comfortable.


-Comfortable walking support

The high waist with extra padding combined with the padded groin straps ensure a high degree of comfort and safety.



Molift Ambulating Vest. Part numbers

		waist dim.
S	Part No. 3036100	23" - 29"
M	Part No. 3036200	29" - 39"
L	Part No. 3036300	39" - 58"

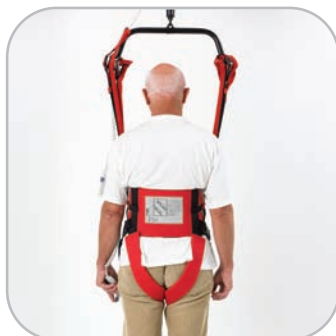

Twin belt system and padded groin straps.

Self-tightening system.

Molift Easy Raiser

Molift Easy Raiser sling is padded and designed for active lifts in combination with Molift Quick Raiser. Providing a perfect dorsal support, it is the ultimate sling for transfers from sitting to sitting, stand-up training and toileting situations.

-Quick and easy

The sling is easily fitted and has a "marine" rope fixing. The design of the sling makes it perfect for toileting.

-Safe and comfortable

Molift Easy Raiser has an anti-slip surface on the inside to prevent the sling from sliding up. It has a high back to distribute pressure evenly, and additional padding under the arms. The waist belt gives extra waist support and stability.



Molift Easy Raiser. Part numbers

		waist dim.
XS	Part No. 3020003	21" - 23"
S	Part No. 3020000	23" - 29"
M	Part No. 3020001	29" - 39"
L	Part No. 3020002	39" - 58"
XL	Part No. 3020004	58" - 68"

Molift Easy Raiser with support

		waist dim.
XS	Part No. 3019003	21" - 23"
S	Part No. 3019000	23" - 29"
M	Part No. 3019001	29" - 39"
L	Part No. 3019002	39" - 58"
XL	Part No. 3019004	58" - 68"

Separate support for Molift Easy Raiser Support

XS/S/M	Part No. 3019011
L/XL	Part No. 3019012



Perfect for stand-up training.



High back for optimal weight distribution.



Molift Easy Raiser Support when extra support is necessary.




Sling size guide

All Molift Slings have different color codes indicating the size. This makes finding the right sling easier. Normal slings for adults range from S to XXL, while pediatric lifting slings are sizes XS and XXS. However, as Molift always focus on comfort and safety, we strongly recommend individual testing and selection.

Size:	Colour guide:	Recom. working load (lb):	Safe working load (lb):
XXS	Pink	26 - 37 lbs	660 lbs
XS	Light blue	37 - 55 lbs	660 lbs
S	Red	55 - 110 lbs	660 lbs
M	Yellow	99 - 209 lbs	660 lbs
L	Green	198 - 352 lbs	660 lbs
XL	Blue	352 - 528 lbs	660 lbs
XXL	White	506 - 660 lbs	660 lbs

Measurements for Molift Easy, Molift Basic and Molift Basic Bathing (inches):

Size:	XXL	XL	L	M	S	XS	XXS
Shoulders (A)	33"	29"	25"	21"	19"	18"	16"
Hips (B)	51"	45"	39"	33"	29"	25"	21"
Back (C)	26"	25"	25"	24"	23"	22"	21"

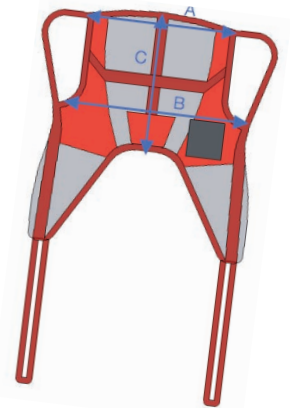

Washing and maintenance

All Molift slings are, unless otherwise indicated on the sling label, washable at 80°C/170°F and autoclaveable at 85°C/185°F for 30 minutes. The plastic insert of head support must be removed before washing. Inspect the slings regularly, especially after laundering. Check carefully for wear and damage to seams, fabric, straps and strap loops. Torn, cut, frayed or broken slings can fail, resulting in serious injury to the user. Always use slings in good condition. Discard and destroy old, unusable slings. When in doubt, please contact the manufacturer or the supplier.

Sling compatibility

Molift allows slings manufactured by companies other than ourselves to be used on our equipment, providing that the slings are of the loop type attachment and a documented risk assessment is carried out by a competent person.

We do, however, recommend our own slings, the main reason for this is, unlike other hoist manufacturers, we have a 4-point suspension (spreader bar) which is designed to be used in conjunction with our slings and not slings designed for a 2-point spreader bar where the clients' position and the stress points of the sling may be compromised.



Accessories and spare parts

Extension loops (2pcs)
Part No. 3013000



Sling bag (folded)
Part No. 3048000



Sling bag (open)
Part No. 3048000



Cushion (for Basic Sling and old Easy Sling w/head support)
Part No. 3006999



Plastic insert
for slings with head support
Part No. 3022601-3022707



Separate support
for Easy Raiser Support
Part No. 3019011-3019012



Important!

Lifting and transferring a patient always involves a degree of risk. Study the Operator Manual and instructions for the accessories carefully. A complete understanding of the contents of the instructions is essential. Only trained personnel should use the equipment. Make certain that the lift accessories are for the lifter you are using and for that particular lifting situation. Exercise caution and care when using the lifting equipment and the accessories. As a caregiver, you are responsible for the safety of the patient. You must therefore be well informed about the patient's condition and ability to manage the lifting situation.